



IN THIS ISSUE

- Printacall Home Room Loop System
- Sennheiser™
- Amplified phones
- Phonic Ear™ personal communicator
- Phonic Ear™ Home Hearing Aid System
- Nokia Loopset for mobile phones
- Specialised amplified sound system for the classroom
- TTY communication options from Telstra

“Improving the life of Deaf and hearing impaired people since 1984.”

Improve your TV Listening

A Printacall Home Room Loop System brings the sound to you


A home room loop enables a person wearing a hearing aid with a T switch to access the sound from a TV, radio or other equipment via the loop amplifier. The loop minimises distracting background noise from within the room and from outside.

The volume of sound from the TV or radio through the room loop to the user is at a more natural level. The whole family can enjoy the sounds of entertainment without uncomfortable noise levels.

The user of the room loop is free to move anywhere within the looped room and enjoy the sounds at their individual volume level, without interfering with other people.

PRINTACALL LOOP SYSTEMS

Does your church or meeting hall have an audio loop system? Printacall is a specialist in the design and installation of powered audio induction loop systems with over 20 years of experience. Our skilled technicians have installed hundreds of loop systems into public access buildings.

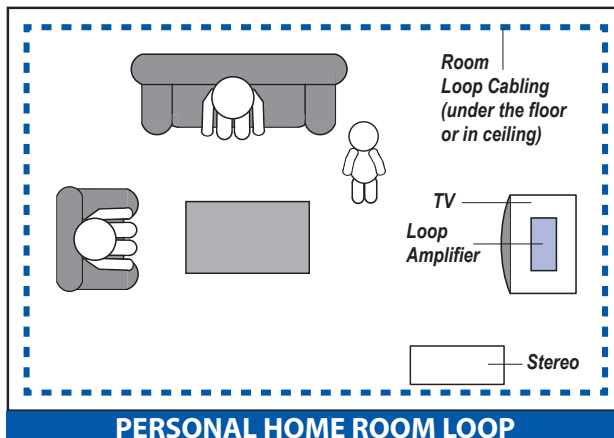

Set 810 Mono/ Stereo Infra-Red System

Turn up the volume when watching TV without disturbing others. Connect the transmitter to the audio output socket on your TV,

then adjust the volume on the headset.



Also available in a neck loop style to deliver high-level quality sound from the TV or radio to your hearing aid with a "T" switch.

* Trademark of Sennheiser electronic GmbH & Co KG.

Products

FEATURED THIS MONTH Oricom TP100*

With extra volume and speech control, the Oricom TP100* is ideal for the hard of hearing. It also boasts a flashing ring indicator so that users can see the phone is ringing.

FEATURES

- Receiving volume control
- Receiving tone control
- Hearing aid compatible
- Sending volume control
- Flashing ring indicator
- Big button keypad
- Easy grip handset
- Ringer volume control
- 3 special one touch memories & 9 one touch memories



* Trademark of Oricom International Pty Ltd.

What is the telecoil in the hearing aid used for?

A hearing aid with an in-built 'telecoil' or 'T switch' affords the wearer access to a range of benefits derived from assistive listening devices - equipment specifically designed and manufactured for compatibility with a hearing aid T switch.



In a conference situation, or in a personal conversation with one person or within a group, the T switch may be coupled with an audio loop assisting you to hear without the interference of background noise.

Nokia Loopset LPS-4

Use your Nokia digital mobile phone without interference



FEATURES

- Used with hearing aid with "T" switch
- Receives power from the phone, so it doesn't require batteries
- Integrated Answer/End button for convenient operation
- Automatic gain control to ensure an optimal audio level
- Nokia model compatibility data available

* Trademark of Nokia Corporation.

Phonic Ear* Easy Listening System

Brings amplified sound direct to your ears!

Transmitter



The person who is talking wears the lapel microphone and transmitter.

A wireless and completely portable system with individual volume control for use indoors or outdoors.

Can be used in restaurants, in cars and in meeting rooms.

Receiver

Also available in Earbuds, Headphones or Neck Loop (for hearing aid wearers with T-switch) styles.



* Trademark of Phonic Ear Inc of USA. Phonic Ear / Logia Denmark

The person who needs assistance with their hearing uses the receiver. The receiver picks up the transmitted signal sent by the transmitter.



Have you received this newsletter in the post and would prefer to receive it via e-mail?

Simply send us an email to let us know: sales@printacall.com.au




Products

Be alerted to the door and phone using a system which you can add to as your needs change

Ask us about other accessories to alert you to various alarms or your baby crying to complete the system.



Place the "Detect" near the sound you want to be alerted to.

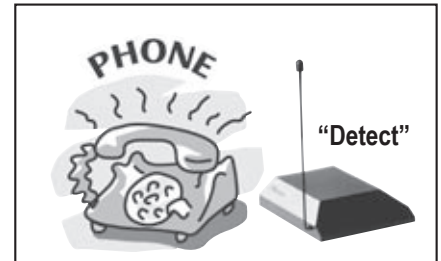

The "Detect" will pick up the sound and send it to an alerting device such as the "Flash" (a light device which

flashes a halogen lamp) or the "Call" (a small loudspeaker which emits a pre-selected audio sound).

(Other types of alerting devices are available)



Phonic Ear* Wake

A bedside alarm and alert device with under pillow vibrator to alert you to the alarm clock or other sounds transmitted from the "Detect".

FEATURES

- Reacts only to frequencies above 600Hz - minimum risk of false alarms
- Microphone detects sound from alarm clock and transforms it into vibrations
- Connects to 240V

* Trademark of Phonic Ear Inc of USA.
Phonic Ear / Loggia Denmark



ORDERING OPTIONS

BY PHONE
02 9809 2392

BY TTY
02 9809 1283

BY POST
2 Doig Ave., Denistone
East, NSW, 2112

BY FAX
02 9809 2345

BY EMAIL
sales@printacall.com.au


For more information visit
www.printacall.com.au

Frontrow™ Pro

Frontrow™ Pro sound field system can benefit your students and teachers!

- Improve the learning environment.
- Engage students and improve attentiveness in the classroom.
- Foster control and orderliness, increase on-task behaviour.
- Boost comprehension, improve test scores.
- Reduce vocal strain and fatigue for teachers.
- Support students suffering from Otitis Media (middle ear infection).

A typical sound field system consists of an infrared emitter and a wireless microphone worn by the teacher and 4 speakers located around the 4 walls of the room which actually make the "Field" integrated with specially designed sensors. Amplifying the teacher's voice raises the signal to noise ratio which improves speech intelligibility and reduces vocal strain. An extra hand held microphone can be provided for students to pass around, so they can be amplified too.

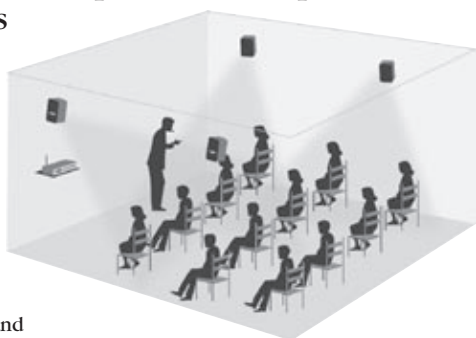

For further information contact Rebecca or Jodie who will be pleased to answer your questions and arrange for you to receive an information pack.

www.soundfield.com.au

Your local distributor

TTY communication options from Telstra

Teletypewriters (TTYs) allow customers to send and receive text conversation over the telephone network. TTYs provide the flexibility to conduct text conversations with other TTY users, or to use the National Relay Service to conduct text to voice calls to other callers.

Telstra offers eligible customers the option of renting one of three TTY models under its extensive Disability Equipment Program.

For further information please contact:

Telstra's Disability Enquiry Hotline
FREECALL* 1800 808 981 (TTY),
FREECALL* 1800 068 424 (Voice),
or Email: Disability@online.telstra.com.au

*A FREECALL from most fixed phones. Printacall is contracted to Telstra to supply, distribute and maintain TTYs.

Uniphone 1150*

The telephone for everyone, combination TTY, telephone and amplified phone all in one!



*Trademark of Ultratec Inc., Madison, WI USA

FEATURES:

- Turbo Code™ with interrupt facility
- Voice carry-over (VCO) and Hearing carry-over (HCO)
- Built-in ring flasher for incoming calls
- Hearing-Aid coupler
- Volume control on handset
- Voice announcer. ASCII
- 48 character LCD display

PRODUCT SPECIFICATIONS AND AVAILABILITY MAY CHANGE WITHOUT PRIOR NOTICE



PH: 02 9809 2392

- TTY: 02 9809 1283
- EMAIL: sales@printacall.com.au
- FAX: 02 9809 2345
- WEBSITE: www.printacall.com.au
- 2 Doig Ave.,
- Denistone East,
- NSW, 2112